

Surgical Technique

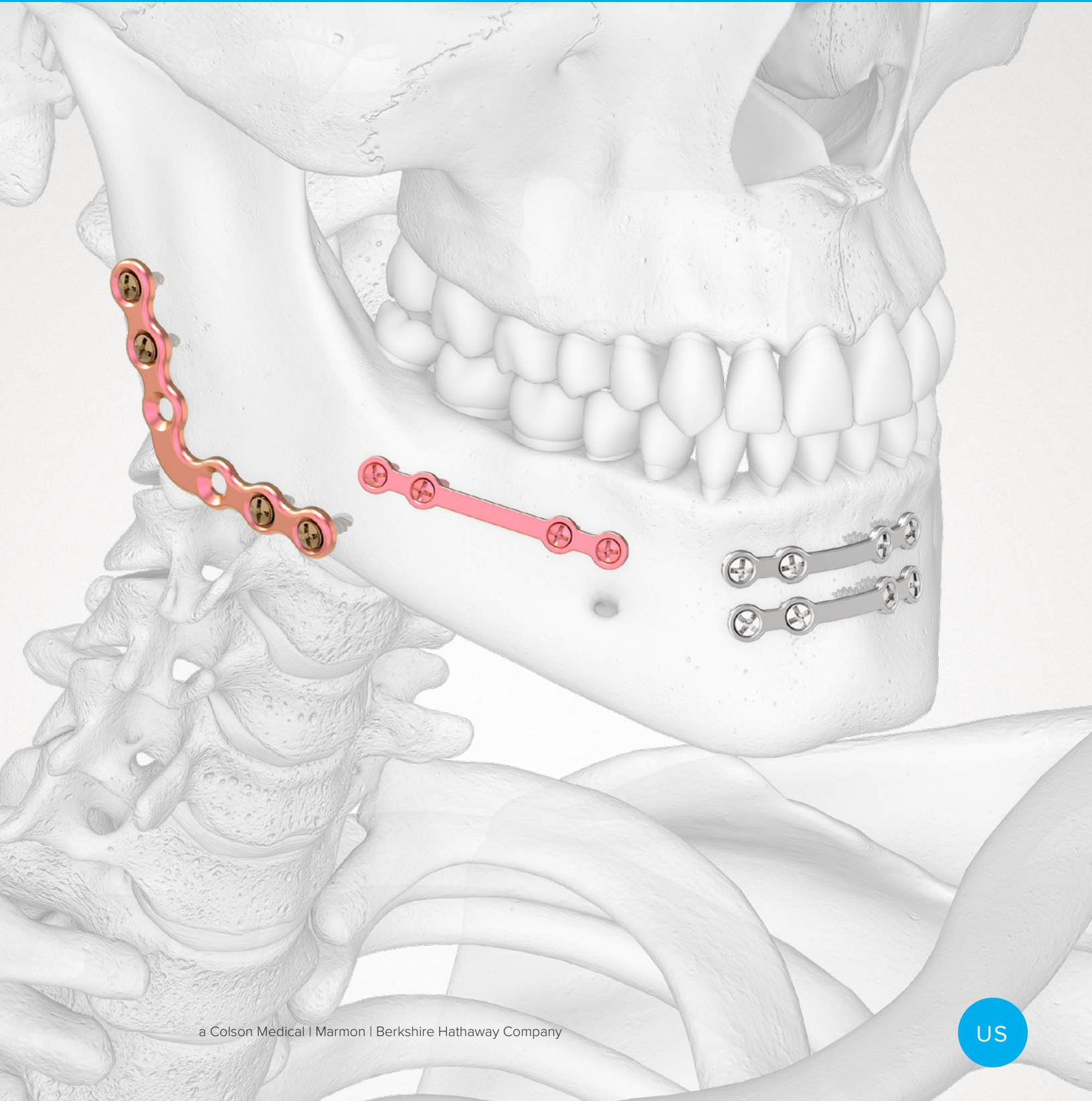




Table of Contents

■ MFX™ Introduction.....	2
■ System Overview	
• Blocks.....	3-4
• Implants.....	5-8
• Instruments.....	9-11
■ Indications.....	12
■ General Surgical Technique Guide	
• Pre-operative planning.....	13
• Fracture reduction.....	14
• Select implants.....	14
• Templates.....	15
• Cutting option.....	15–16
• Contouring options.....	17–20
• Plate placement.....	20
• Pilot hole drilling.....	21–22
• Drilling options.....	22
• Determining screw length.....	23
• Screw insertion.....	24–25
• Placement of remaining screws.....	25–26
• Close procedure.....	26
■ Specific Techniques	
• Transbuccal technique.....	27–31
• Locking technique.....	32–34
• Angulated locking technique.....	35–38
• Reduction plates and forceps.....	39–43
• Bone resection technique.....	44–45
• Lag screw technique.....	46–48
• MMF technique (Maxillomandibular Fixation)...	49–53

Introduction


OsteoMed introduces the MFX™ Rigid Fixation System as the next step in the evolution of mandibular rigid fixation and reconstruction.

The system applies the best of OsteoMed innovation and design to a proven understanding of established, traditional mandible fixation & reconstruction techniques and principles.

The MFX™ system is intended for fracture fixation during orthognathic reconstruction, mandibular reconstruction and surgery involving osteotomies and trauma.

The system is simple, versatile, and dynamic in its application, incorporating the following screw technologies: standard, Auto-Drive®, locking, angulated locking, angulated locking Auto-Drive®, safety, MMF and lag screws.

The implants are made from titanium alloy (ASTM F-136) and/or commercially pure titanium (ASTM F-67) which are MRI compatible.

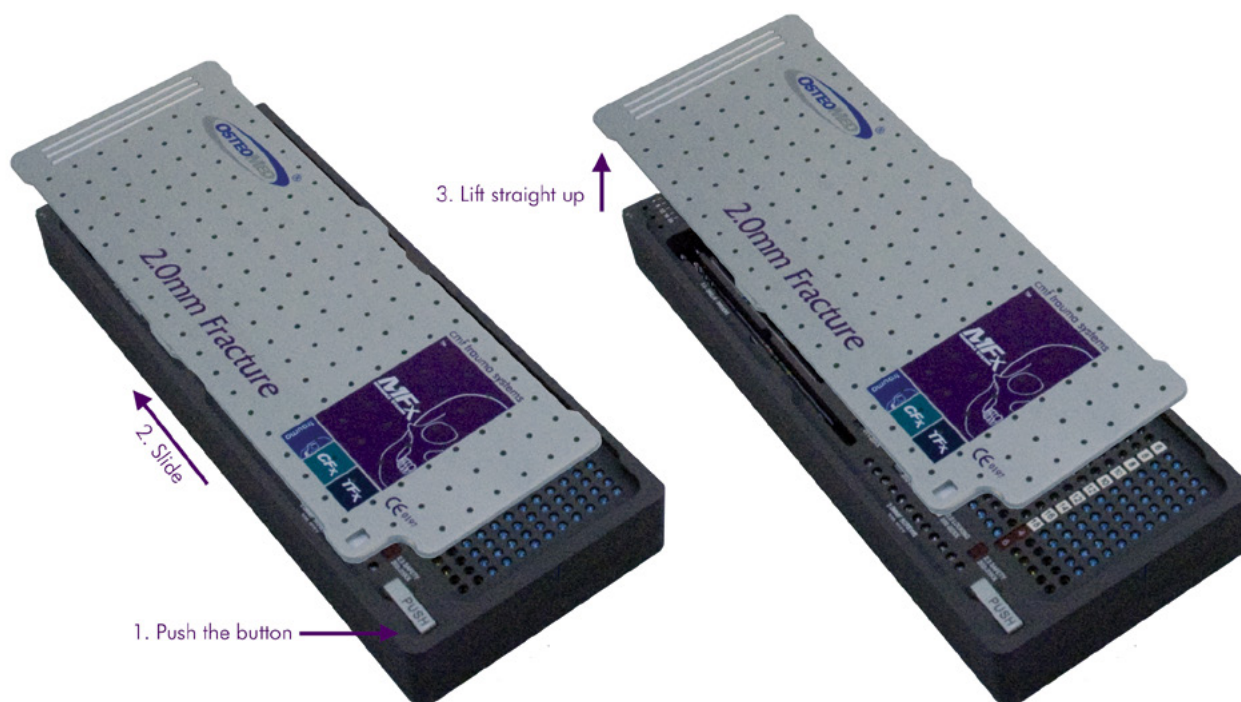

Note: Remove lid from block by holding down the button, pull and lift to open



System Overview


Blocks

Standard Blocks


- MFX™ 2.0 mm Fracture Organizer Block 220-0721



- MFX™ 2.4 mm Fracture Organizer Block 220-0722

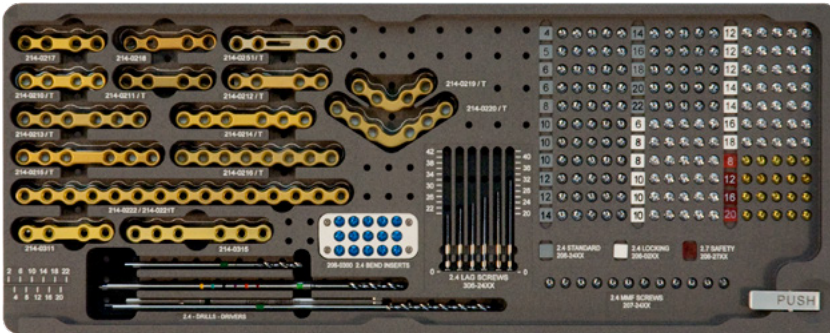

- MFX™ 2.4 mm Reconstruction Organizer Block 220-0723


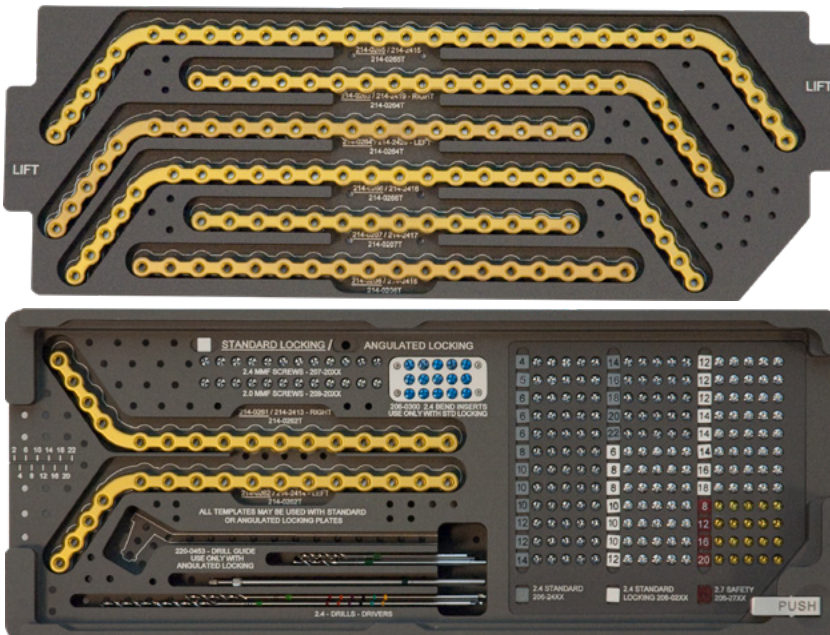


System Overview



Blocks

Optional Blocks

- MFX™ 2.0/2.4 mm Façture Angulated Locking Organizer Block 220-0631

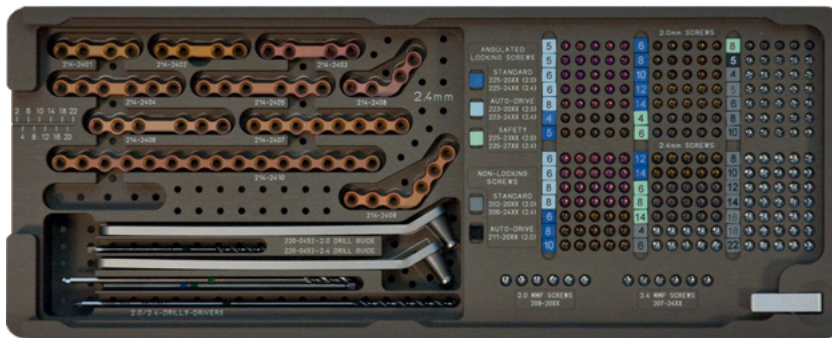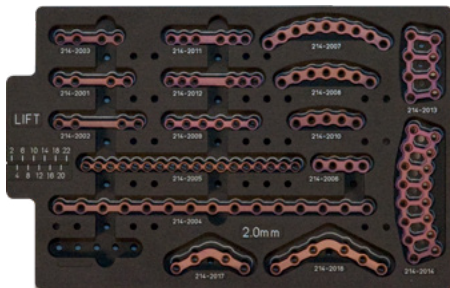

- MFX™ 2.4 mm Reconstruction Angulated Locking Organizer Block 220-0632

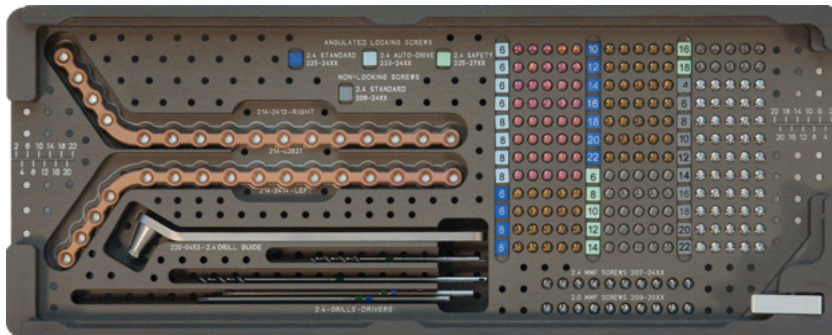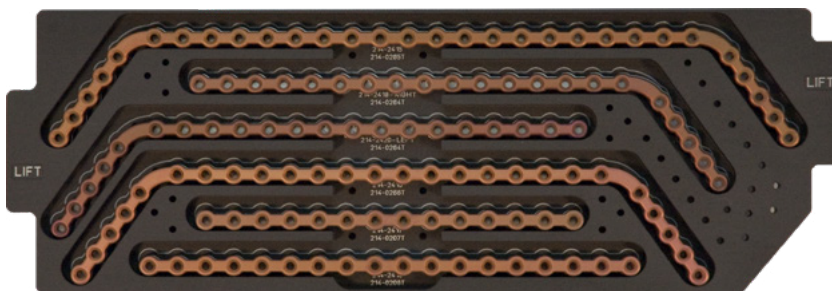

System Overview

Screws*

Standard Screw

- 2.0 mm 202-20XX (4–22 mm)
- 2.4 mm 206-24XX (4–22 mm)



Standard Auto-Drive® Screw

- 2.0 mm 211-20XX (4–8 mm)



MMF Screw

- 2.0 mm Auto-Drive® 209-20XX (8–14 mm)
- 2.4 mm 207-24XX (12–18 mm)



Standard Safety Screw

- 2.3 mm 202-23XX (4–6 mm)
- 2.7 mm 206-27XX (8–20 mm)



Standard Locking Screw

- 2.0 mm 202-02XX (6–8 mm)
- 2.4 mm 206-02XX (6–18 mm)



Lag Screw

- 2.4 mm 306-24XX (20–38 mm)



Angulated Locking Auto-Drive® Screw

- 2.0 mm 223-20XX (5–8 mm)
- 2.4 mm 223-24XX (6 mm and 8 mm)



Angulated Locking Standard Screw

- 2.0 mm 225-20XX (4–18 mm)
- 2.4 mm 225-24XX (6–22 mm)



Angulated Locking Safety Screw


- 2.3 mm 225-23XX (4–18 mm)
- 2.7 mm 225-27XX (6–22 mm)




System Overview

Plates*

- 2.0 mm Fracture Block

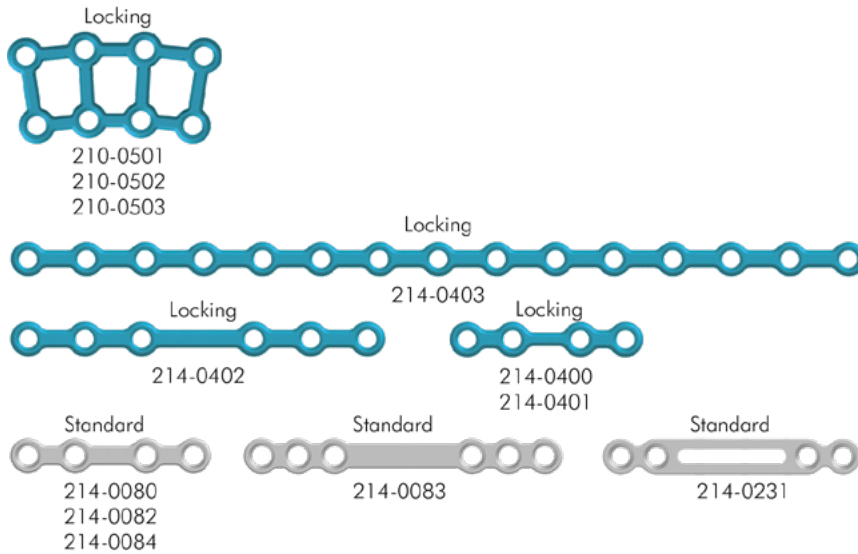

- 2.4 mm Fracture Block

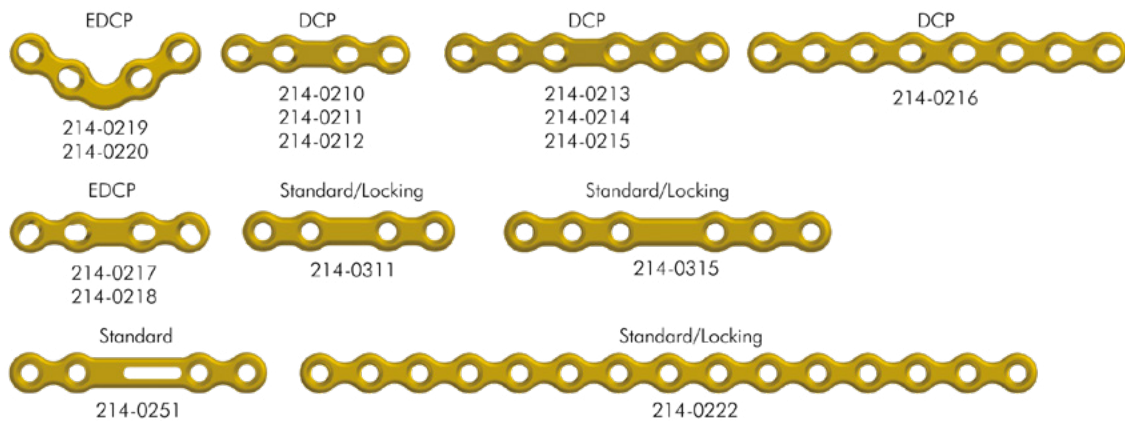





System Overview

Plates*

- 2.4 mm Reconstruction Block

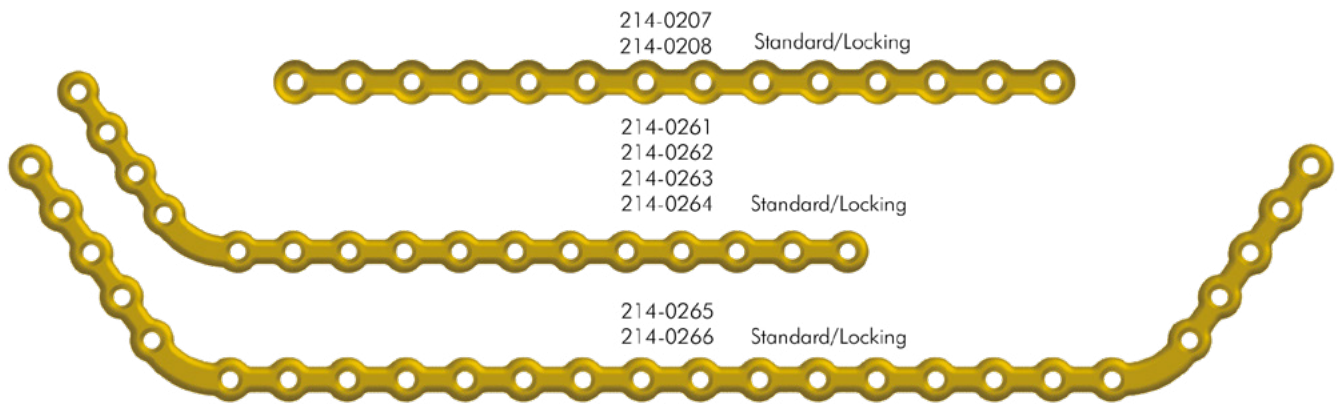


- 2.0/2.4 mm Fracture Angulated Locking Fixation System

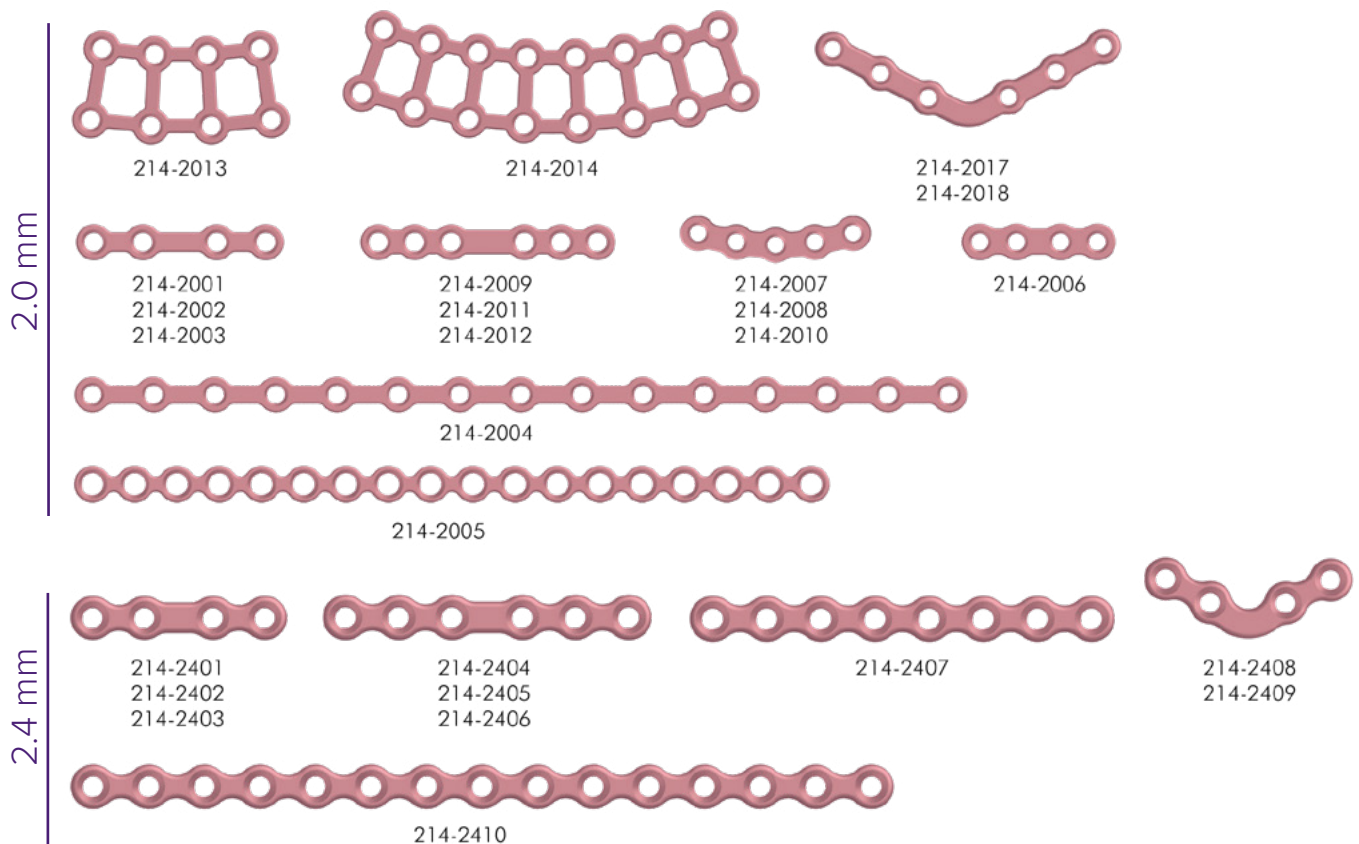

System Overview

Plates*

- 2.4 mm Reconstruction Angulated Locking Fixation System

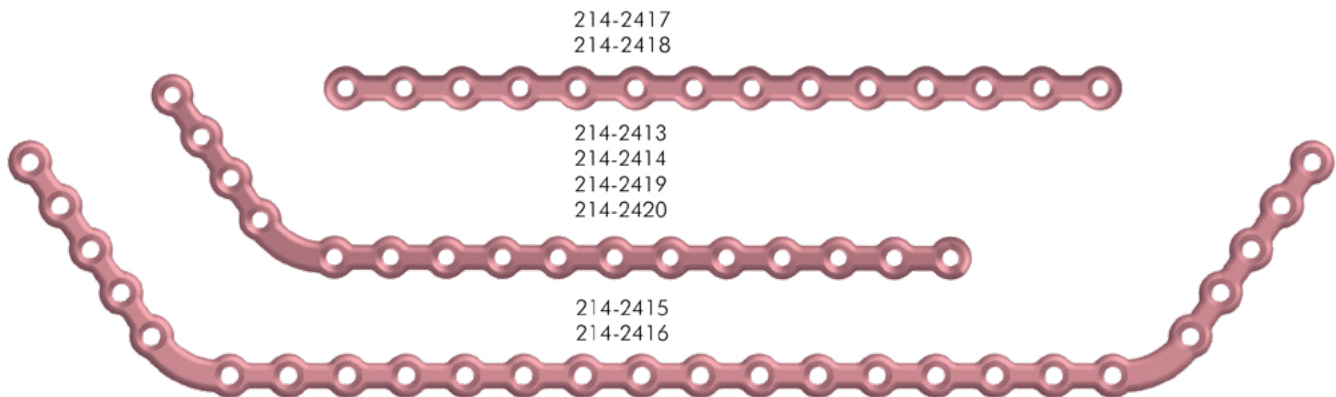


Screws Tags

- 220-0616-XX (Black): Tag, Auto-Drive® Screw
- 220-0617-XX (Bone): Tag, Locking Screw
- 220-0614-XX (Gray): Tag, Standard Screw
- 220-0615-XX (Orange): Tag, Safety Screws
- 220-0613-XX (Dark blue) MFx™ Tag, Standard Angulated Locking Screw
- 220-0618-XX (Light blue) MFx™ Tag, Auto-Drive® Angulated Locking Screw
- 220-0619-XX (Light green) MFx™ Tag, Angulated Locking Safety Screw



System Overview

Instruments

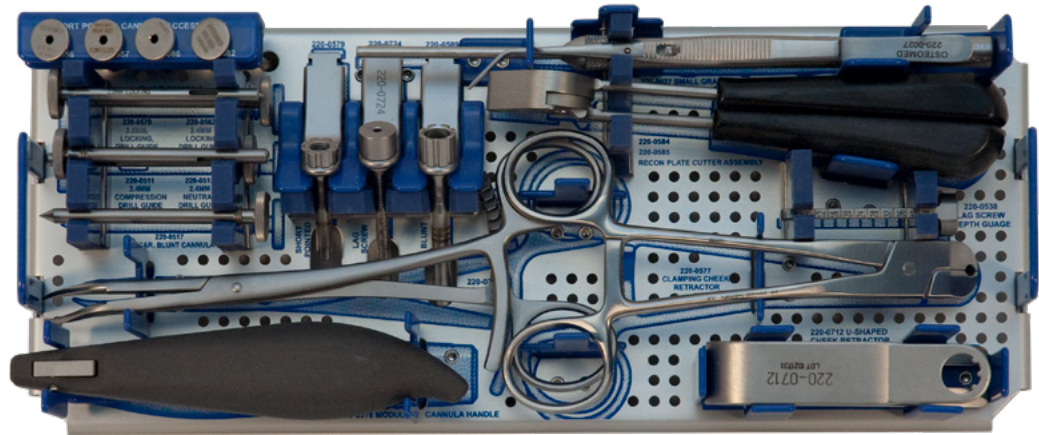

Tray 1



220-0578 Modular Handle, Cannula



220-0712 Cheek Retractor, U-Shaped



220-0538 Sliding Length Gauge, Lag Screw Depth Gauge



220-0579 Modular Cannula, Short Pointed



220-0589 Modular Cannula, Mandible/Blunt




220-0231 Reduction Forceps



220-0577 Cheek Retractor, Clamp

220-0724 MFX™ Lag Screw Cannula

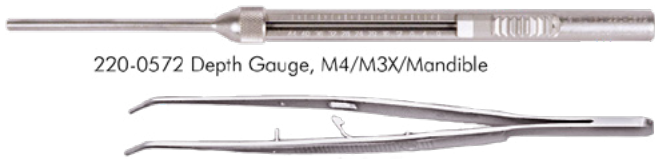

220-0572 Depth Gauge, M4/M3X/Mandible



220-0584 Plate Cutter, Reconstruction, Unflanged



220-0586 2.4mm Neutral Drill Guide



220-0556 2.4mm Locking Drill Guide CMF



220-0557 2.0mm Locking Drill Guide CMF



220-0582 Trocar, Short Pointed



220-0511 Drill Guide, Compression, Blunt Cannula



220-0513 Drill Guide, Neutral, Blunt Cannula



220-0562 Drill Guide, 2.4mm Locking Plate, Blunt Cannula



220-0570 Drill Guide, 2.0mm Locking Plate, Blunt Cannula




220-0517 Trocar, Blunt Cannula



220-0585 Plate Cutter, Reconstruction, with Flange

System Overview

Instruments

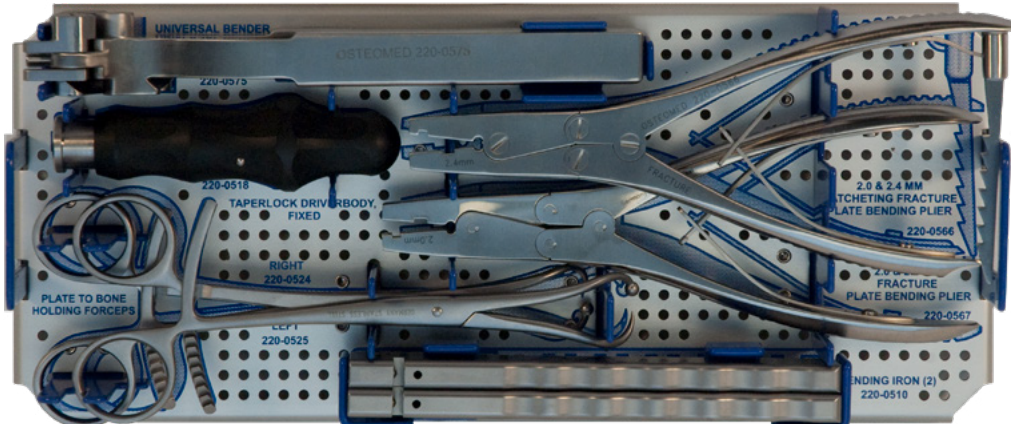

Tray 2



220-0525 Plate to Bone Holding Forceps, Left



220-0524 Plate to Bone Holding Forceps, Right



220-0518 Taperlock Body, Mand, Fixed



220-0575 Universal Plate Bender

(2)



220-0510 Bending Iron, Recon Plates




220-0566
2.0/2.4mm Fracture Plate Bending Plier, Ratchet



220-0567
2.0/2.4mm Fracture Plate Bending Plier

System Overview


Instruments

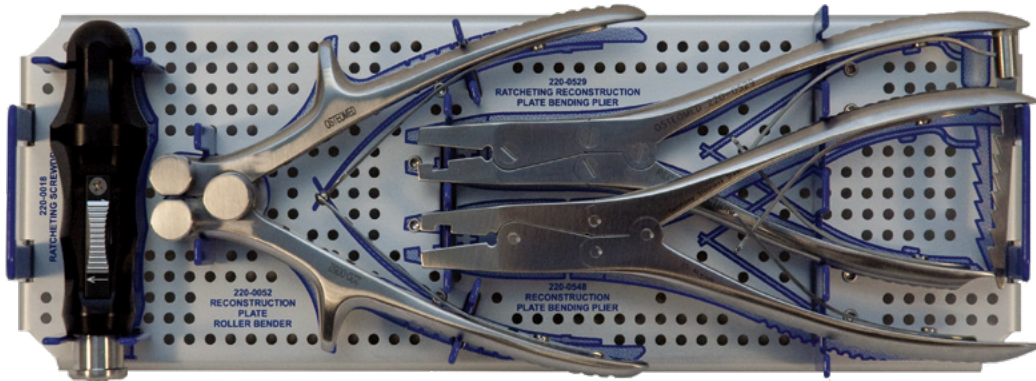

Tray 3



220-0018 Ratcheting Screwdriver

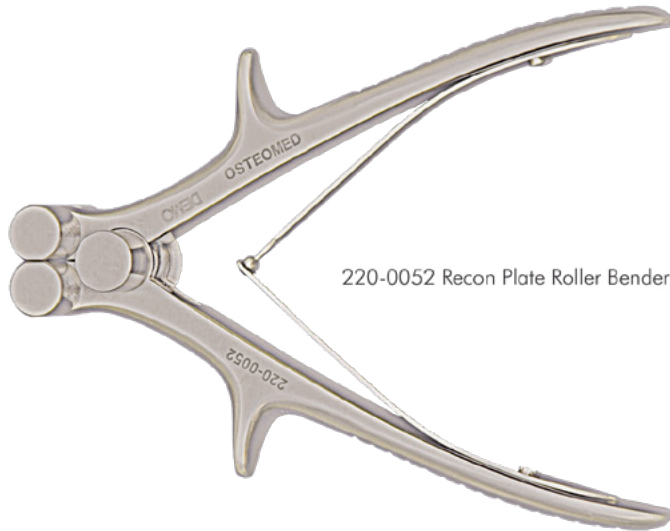

220-0052 Recon Plate Roller Bender



220-0548
Reconstruction Plate Bending Plier, No Ratchet

220-0529
Reconstruction Plate Bending Plier, w/Ratchet

Indications

- The OsteoMed MFx™ Mandibular Fracture / Reconstruction System is indicated for fracture fixation, mandibular reconstruction and surgery involving osteotomies and trauma.
- The OsteoMed 2.0 Locking Plate System is indicated for oral, maxillofacial surgery, trauma, reconstructive surgery and orthognathic surgery (surgical correction of dentofacial deformities).
- The OsteoMed MMF Screws are indicated for temporary ligature and wire lock fixation for temporary constriction and stabilization of fractured bone segments in the oral cavity in conjunction with primary fixation devices.
- The OsteoMed Angulated Locking Fixation System is indicated for mandibular trauma reconstruction, mandibular reconstruction and orthognathic reconstruction.
- The OsteoMed Reduction Plates and Forceps are intended for mandibular body fractures, symphysis fractures and parasyphysis for tension, compression or both.
- The OsteoMed implants, templates and drills are intended for single use only.

General Surgical Technique Guide

Pre Operative Planning

Fractures of the mandible will typically reveal a malocclusion (inability to bite down), pain at fracture site, significant internal bruising, or laceration with bleeding between teeth at the fracture site.

Reduction and stabilization of the mandible fracture is key to successful treatment. The method of management may vary based on the severity, location of the fracture, presence or absence of teeth and patient demographic. Mandible fractures may be treated by closed reduction with wiring of the teeth or open reduction with internal rigid fixation using plates and screws. The technique of closed reduction may involve wiring the teeth for up to 4 to 6 weeks.

Internal rigid fixation when performing open reduction requires exposure of the fracture sites and stabilization with screws and/or plates. Accurate reduction with good stabilization can frequently avoid complications and help to restore the patient's primary occlusion and facial appearance. In the following pages we will illustrate open reduction techniques¹.

General Surgical Technique Guide


1. Reduce the fracture.




Use Reduction Forceps (220-0231):

- Pre-drill anchor holes within fracture segments and insert Reduction Forceps.
- Squeeze Reduction Forceps until fracture is reduced.

2. Select plates.

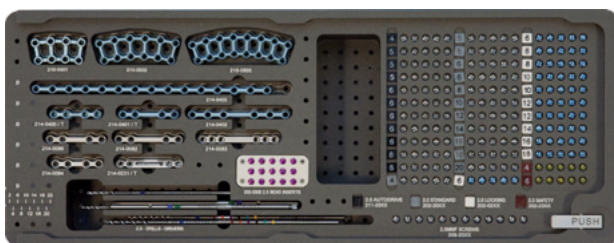

MFx™ 2.0 mm Fracture Organizer Block
(220-0721)

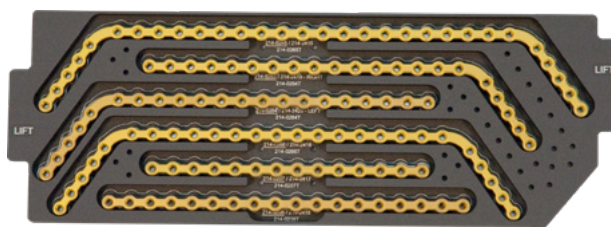

MFx™ 2.4 mm Reconstruction Organizer Block
(220-0723)

Technique Tip: Ensure the implants achieve the following:

- Adequately bridge the defect
- Provide stable fixation
- Fit structure of anatomical region
- Return form and function to fractured bone

General Surgical Technique Guide

3. Select the template if applicable to plate selection.



Contour the template to match the defect.



4. Cut the plate if necessary.




(220-0584)



(220-0585)


Use Round Reconstruction Plate Cutter Assembly. (220-0584, 220-0585):

- Place plate into handle. (220-0585)

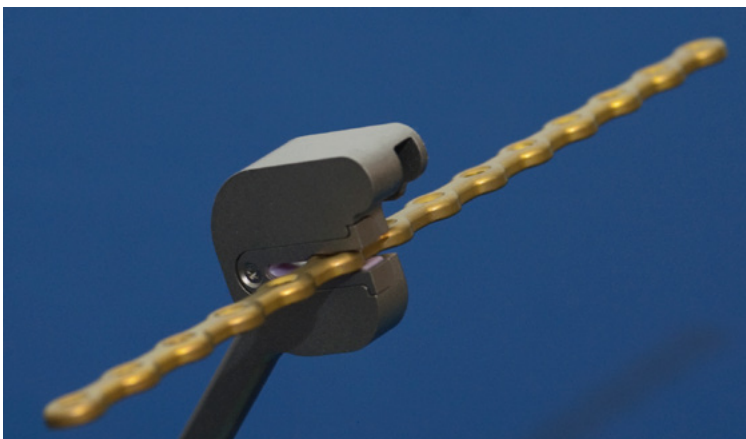

General Surgical Technique Guide

- Slide second handle into first handle until plate is secure.

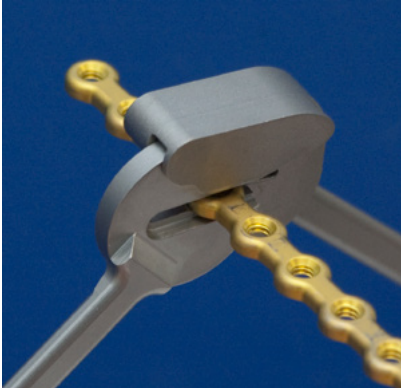

- Rotate handles in towards each other and then back apart, repeat this until plate is separated into 2 pieces.



- Utilize diamond file on the handle to remove sharp edges on the plate. Located on 220-0584.



General Surgical Technique Guide

5. Contour the plate.

Utilizing the contoured template, contour the plate to match the template.

Instrumentation for Fracture Plates:



■ Fracture Plate Bending Pliers (220-0566, 220-0567):

- Use with appropriate plates- 2.0 mm/2.4 mm fracture plates.
- Place plate in Fracture Plate Bending Pliers so the plate fits into recessed area.
- Manipulate plate to desired shape.

Note: All standard locking plates require bending inserts.

Note: Medial holes will not fit into recess.


Note: Anterior flat in bending pliers is to be used to contour medial portion of plates.

Note: If a Reconstruction Plate is selected for highly comminuted fractures or ablative procedures the following instrumentation options are available

■ Reconstruction Plate Bending Pliers (220-0529, 220-0548):

- Select appropriate plates - 2.4 mm reconstruction plates
- Place plate in Reconstruction Plate Bending Pliers so the plate fits in recessed area
- Manipulate plate to desired shape

Note: When using a pre-bent reconstruction plate, screw holes in bend will not recess into pliers

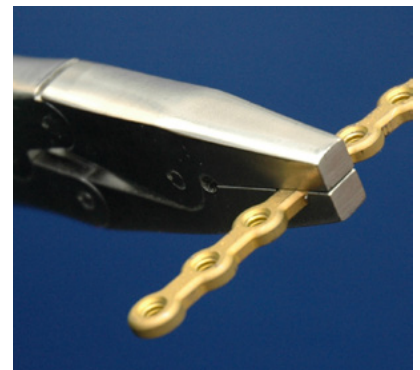

General Surgical Technique Guide

■ Universal Plate Bender (220-0575):



- Place plate in 1 of 3 options on Universal Plate Bender:

Note: *Do not bend directly over the plate hole*

Option 1:


- Place plate between prongs
- Squeeze Universal Plate Bender until plate bends to desired angle



Option 1

Option 2:

- Place plate between full circle and 2 semi-circles
- Squeeze Universal Plate Bender until plate bends to desired angle

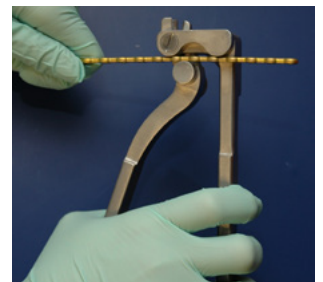

Option 2

Option 3:

- Place plate between semi-circle and concave surface (semi-circle)
- Squeeze until plate bends to desired angle



Option 3

General Surgical Technique Guide


■ Reconstruction Roller Plate Bender (220-0052):

- Select appropriate plates - 2.4 mm reconstruction plates
- Place plate in 1 of 2 options on Roller Plate Bender

Note: *Do not bend directly over the plate hole*

Option 1:

- This option should be used to create bends in a plane perpendicular to the plate holes
- Place plate between 3 cylinders
- Squeeze Roller Plate Bender until plate bends to desired angle

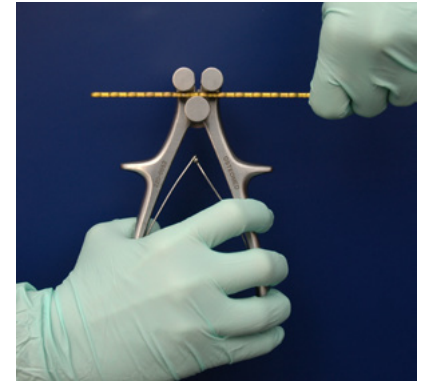


Option 1

Technique Tip: To create smooth continuous bends, bend the plate partially then slide the plate slightly and bend more. Continue this process until the desired shape is acquired.

Option 2:

- This option should be used to create bends in the same plane as the plate holes
- Place plate between 2 cylinders and wedge step
- Squeeze Roller Plate Bender until plate bends to desired angle

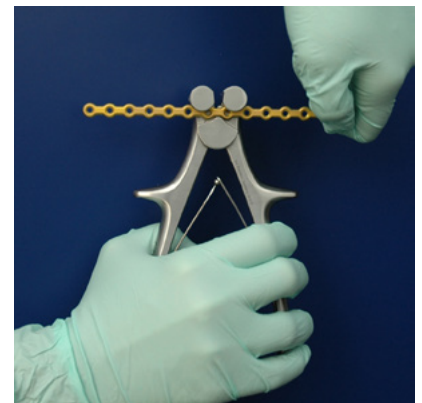

Option 2



Note: *Bends greater than a few degrees should be created by bending at multiple locations*

General Surgical Technique Guide

■ Bending Iron (220-0510):



- Place plate into 1 of 3 recessed areas
- Manipulate plate to desired shape
- Use with appropriate plates - 2.4 mm reconstruction plates

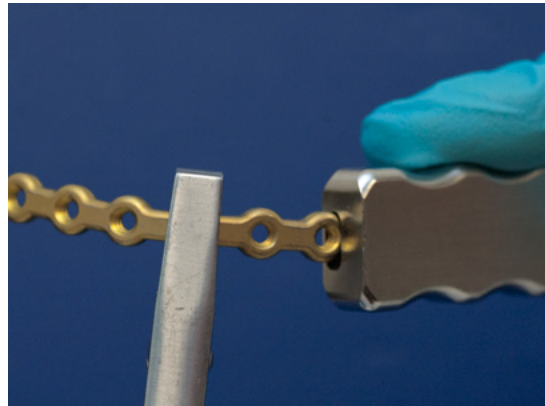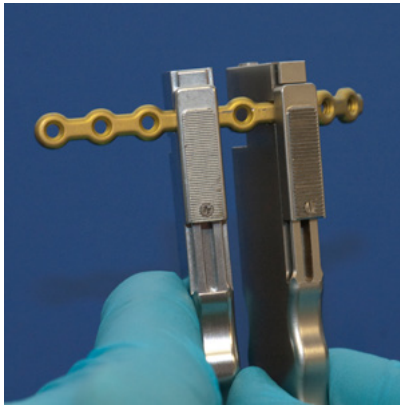

6. Place the plate.



General Surgical Technique Guide

7. Drill Pilot Hole

When a pilot hole is necessary, follow sequence below to achieve desired surgical results.




- Select appropriate diameter drill bit based on screw selection and surgical application.

Note: Blue color bands denote drill for 2.0 mm screws (1.6 mm hole).

Green color bands denotes drill for 2.4 mm screws (2.0 mm hole)


Note: Trocar system requires use of appropriate drill guide.

- Select appropriate drill guide, cannula and attach to Modular Cannula Handle. 220-0578

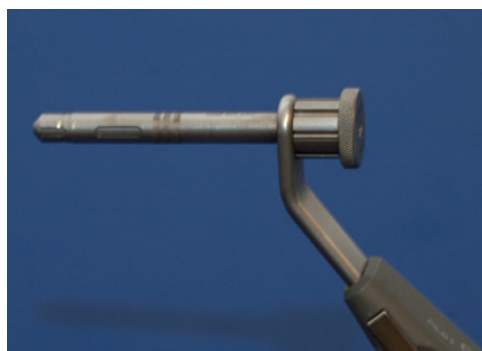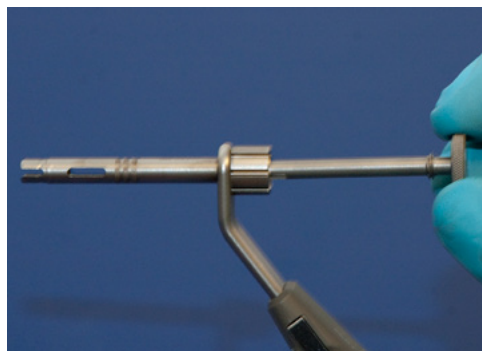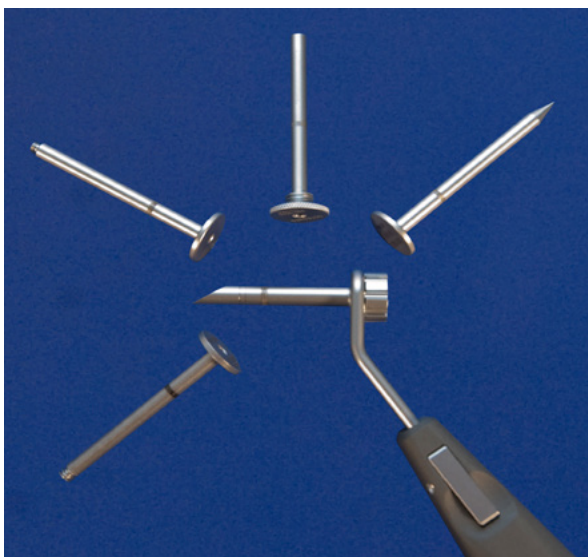

General Surgical Technique Guide

- Drill first pilot holes most proximal to fracture line/osteotomy

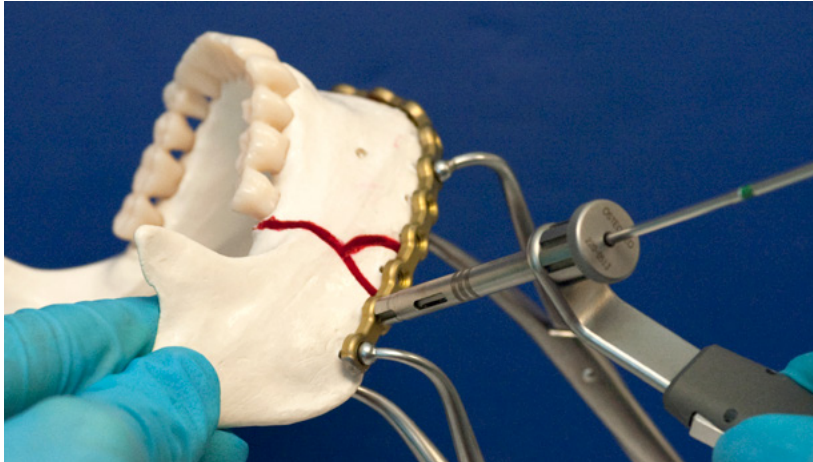

Technique Tip: Rainbow drill legend is located on the back of the organizer block lid.

Technique Tip: Single and triple bands on rainbow drills denote appropriate cannula systems



Note: Pilot hole depth should match or exceed the length of screw. Rainbow drills are in 4 mm increments.

Note: The proximal color bands on Rainbow Drills indicate the depth of pilot holes (220-0251, 220-0252, 220-0255, 220-0256)

Note: Drill speed and torque must follow power system parameters

Note: Use irrigation to prevent bone necrosis

Note: When using a drill guide do not apply a side load on the drill. This may result in friction, which may generate a thermal burn. Axial loading should always be used.


General Surgical Technique Guide

8. Determine screw length

Use Depth Gauge (220-0572) in Mandible:



- Insert Depth Gauge into pilot hole until it engages lingual cortex

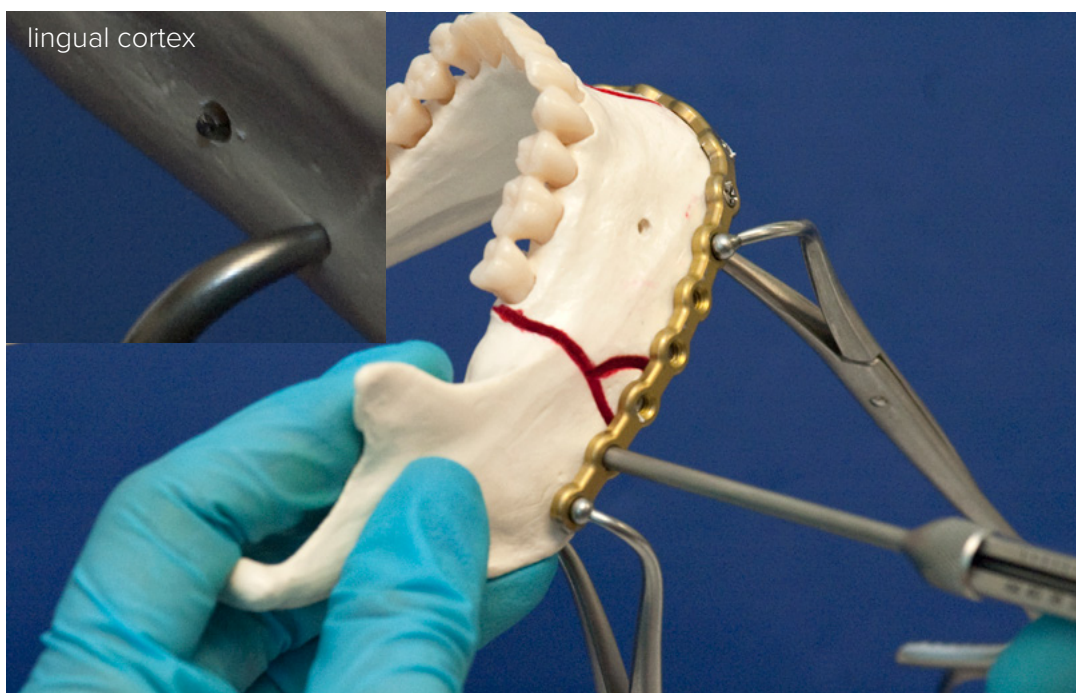

- Determine necessary screw length using scale marked on Depth Gauge.

Note: Always measure screw depth through the plate

Note: When measurement falls between two markings, the longer length should be considered for screw selection

General Surgical Technique Guide

9. Select desired screwdriver body

Screw Driver Options:

■ Ratcheting Screwdriver Body



(220-0018)



■ Taperlock™ Screwdriver Body



(220-0518)

10. Load screw on driver

- Insert appropriate driver shaft into appropriate screw, applying a perpendicular force that engages screw
- Lift straight up to remove screw from Organizer Block

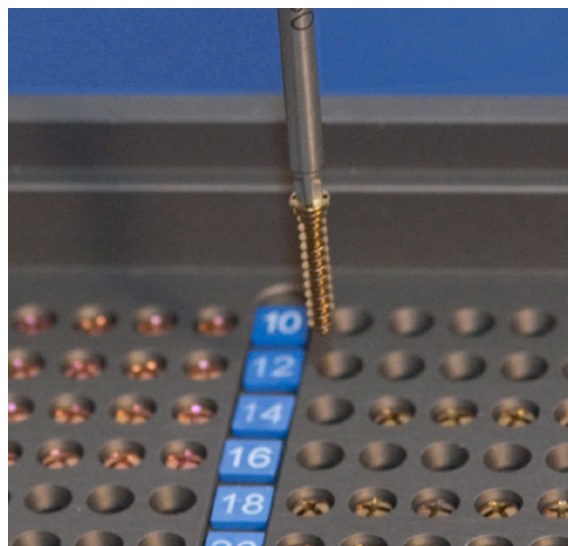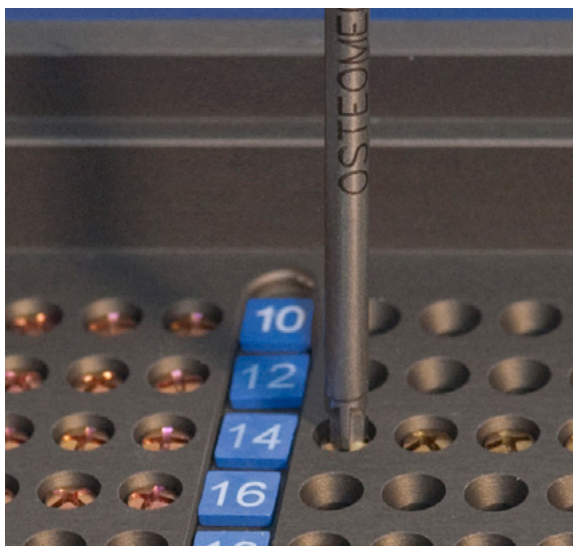


Note: *Organizer Blocks have recessed screw holes to ease loading of screws and are customizable to meet surgeon's screw size needs*

General Surgical Technique Guide

11. Insert screw into bone

Standard:

- Drive first screw into pilot hole most proximal to fracture/osteotomy until screw is flush with surface of bone/plate

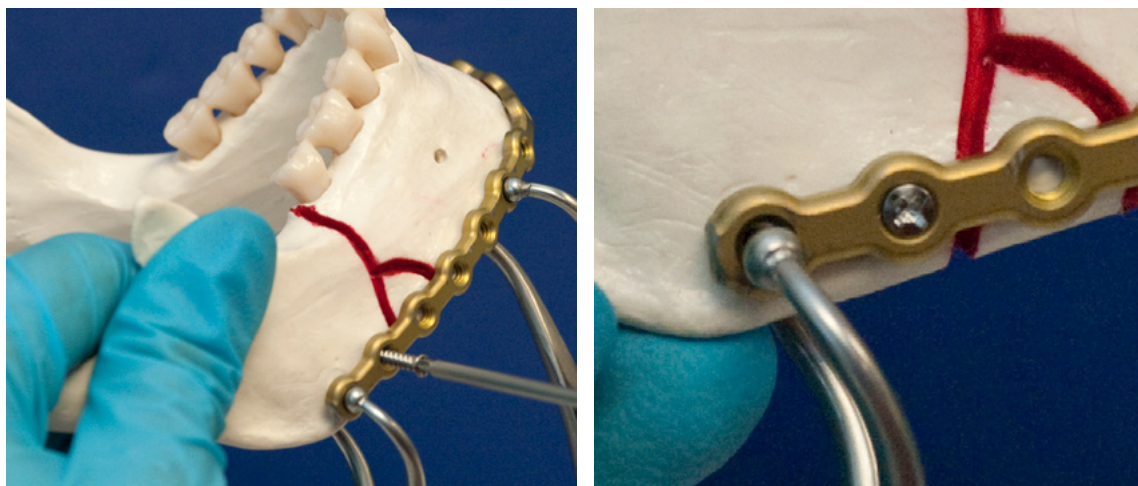

Note: the plate may be held in place using the Plate to Bone Holding Forceps (220-0525, 220-0524)

- Drive second screw into most proximal pilot hole on opposite side of fracture/osteotomy until screw is flush with surface of bone/plate
- Drive remaining screws in same manner in a distal orientation from fracture line/osteotomy



General Surgical Technique Guide

Auto-Drive®:

- Drive first screw into bone, most proximal to fracture.
Note: *Drive the screw perpendicular to the bone until screw is flush with surface of plate/bone*
- Drive second screw into most proximal pilot hole on opposite side of fracture/osteotomy until screw is flush with surface of bone/plate
- Drive remaining screws in same manner in a distal orientation from fracture line/osteotomy

Note: *Higher torque may be required to fully seat screws when using Auto-Drive® screws*

Note: *In high density bone pilot drilling may be necessary*


12. Close per standard practice





Transbuccal Approach Surgical Technique Guide

MFX™ Modular Cannula Drill Guide Selection

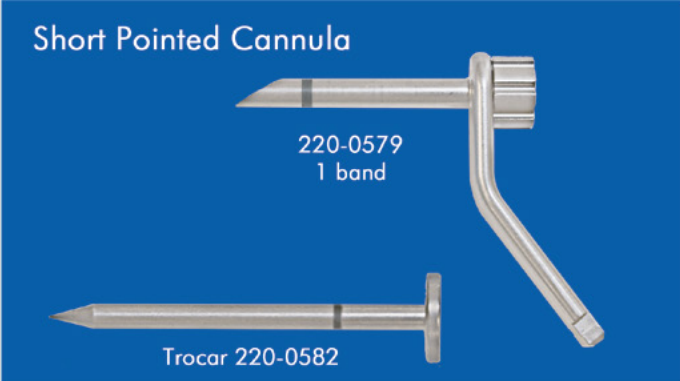

Technology	Drill Guide
Non Locking Plate w/ Standard Screw	220-0586
2.0mm Locking Plate w/ Locking Screw	220-0557
2.0mm Locking Plate w/ Standard Screw	220-0586
2.4mm Locking Plate w/ Locking Screw	220-0556
2.4mm Locking Plate w/ Standard	220-0586
2.0mm Angulated Locking Plate w/ Angulated Locking Screw	220-0586
2.0mm Angulated Locking Plate w/ Standard Screw	220-0586
2.4mm Angulated Locking Plate w/ Angulated Locking Screw	220-0586
2.4mm Angulated Locking Plate w/ Standard Screw	220-0586

Technology	Drill Guide
Non Locking Plate w/ Standard Screw	220-0513
2.0mm Locking Plate w/ Locking Screw	220-0570
2.0mm Locking Plate w/ Standard	220-0513
2.4mm Locking Plate w/ Locking Screw	220-0562
2.4mm Locking Plate w/ Standard Screw	220-0513
2.0mm Angulated Locking Plate w/ Angulated Locking Screw	220-0586
2.0mm Angulated Locking Plate w/ Standard Screw	220-0586
2.4mm Angulated Locking Plate w/ Angulated Locking Screw	220-0586
2.4mm Angulated Locking Plate w/ Standard Screw	220-0586
2.4mm DCP or EDCP plate w/ Standard Screw	220-0511

Note: *Angulated Locking Plates are intended to be used with the Angulated Locking Drill Guide (220-0452 and 220-0453). If they are used with the Modular Cannula System the neutral drill guide must be used.*

1. Select trocar system

Transbuccal Approach Surgical Technique Guide

2. Select cheek retractor

Cheek retractor Options:

- Cheek Retractor, U-Shaped (220-0712)



- Cheek Retractor, Clamp (220-0577)



3. Select appropriate Modular Cannula and attach to Modular Cannula Handle





Cannula Options:

- Short Pointed Modular Cannula (220-0579) - see chart on page 27
- Modular Cannula, Mandible/Blunt (220-0589) - see chart on page 27

Transbuccal Approach Surgical Technique Guide

4. Select appropriate Trocar by matching banding on Modular Cannula with banding on Trocar

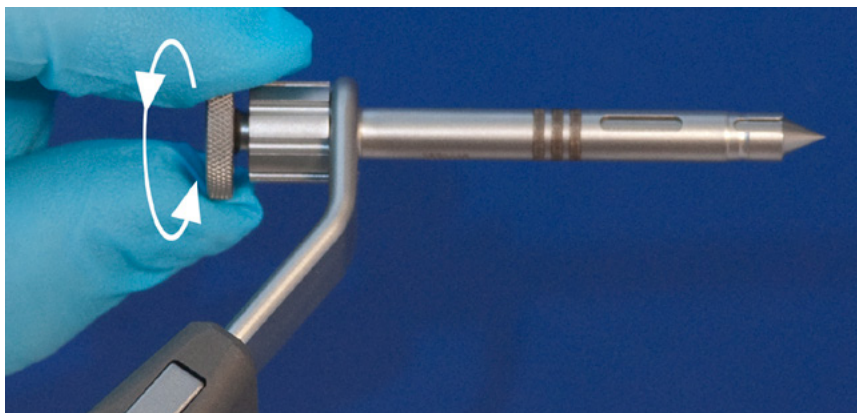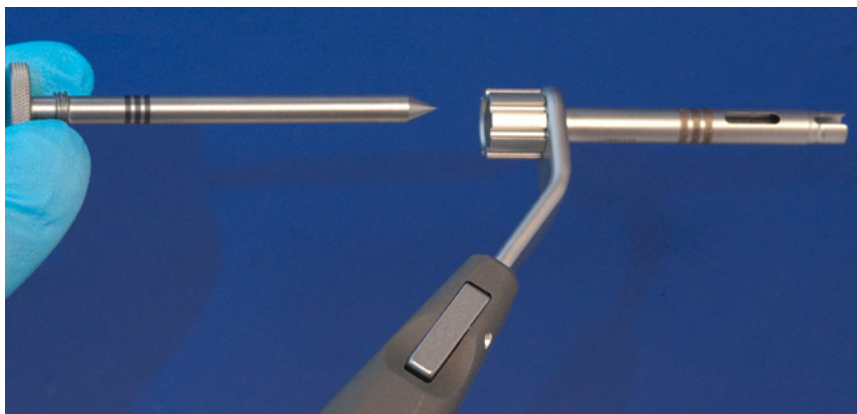

Transbuccal Approach

Surgical Technique Guide

5. Insert Trocar into Modular Cannula and punch through previously made stab incision in cheek

6. Remove Trocar from Modular Cannula

7. Attach Cheek Retractor to Modular Cannula

Note: *Cannula is to never be used as a cheek retractor*



8. Position plate where desired


Transbuccal Approach Surgical Technique Guide



9. Select appropriate Drill Guide by matching banding on Modular Cannula with banding on Drill Guide while selecting appropriate drill guide based on screw/plate technology



10. Insert Drill Guide into Modular Cannula

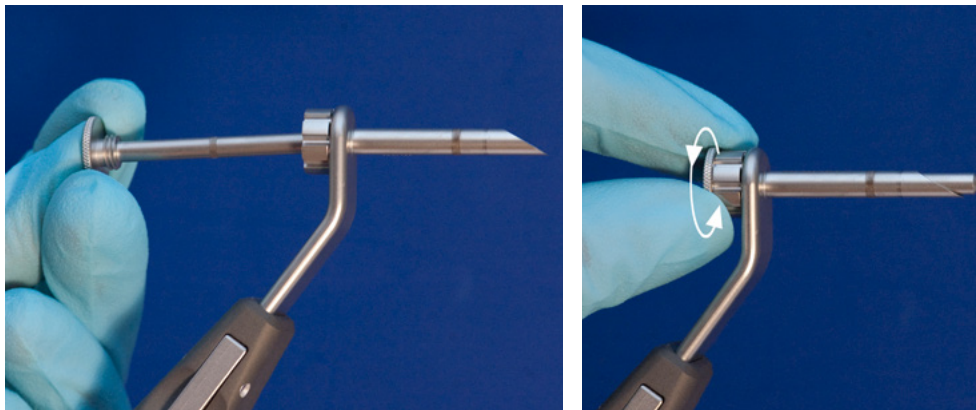


11. Select plate, contour and position on bone
12. Follow instructions for pilot hole drilling
(Reference General STG pg 21–22, Step 7)

Locking Screws

Surgical Technique Guide

1. Expose the defect, select and contour the appropriate template

Note: *If contouring is necessary place Bending Inserts into the fracture plate*

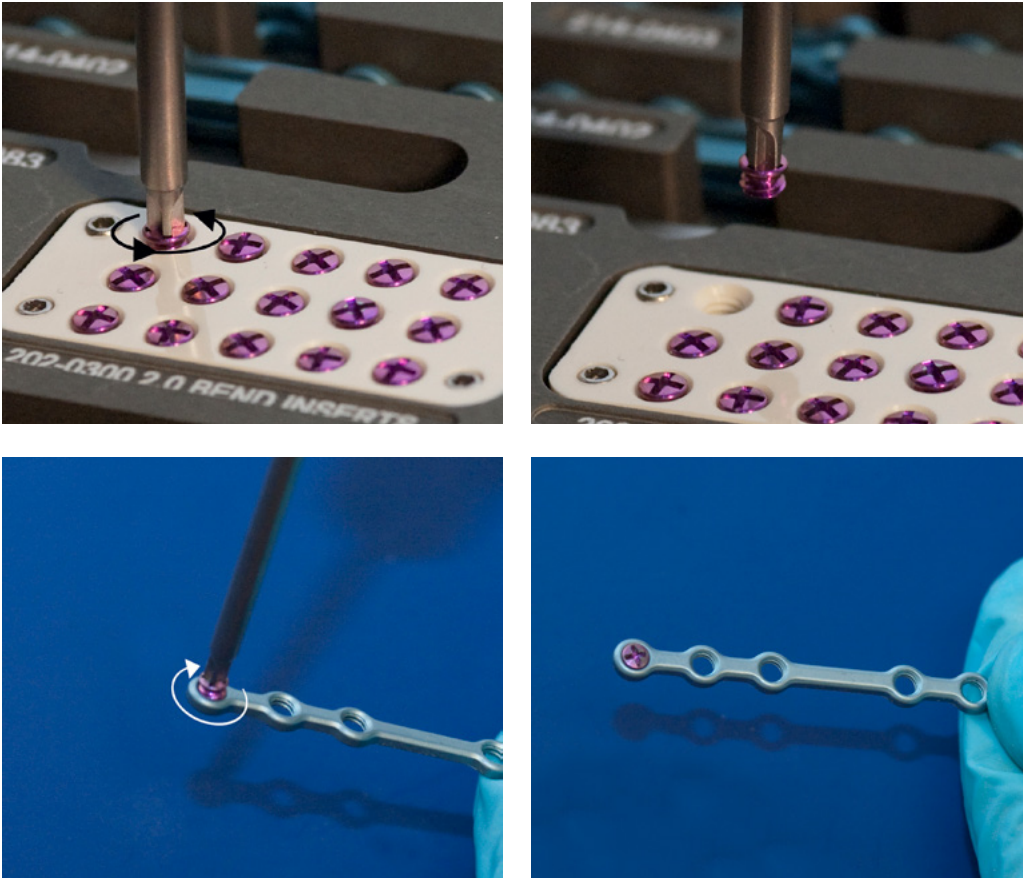

2. Contour the implant using the plate contouring instrumentations to match the template.

(Reference General STG pg 17-20, Step 5)

Note: *remove bending inserts after contouring is completed*

3. Cut the fracture plate, if necessary, using the Round Reconstruction Plate Cutter Assembly (220-0584, 220-0585)

(Reference General STG pg 15-16, Step 4)


Locking Screws

Surgical Technique Guide



4. If transbuccal approach is needed, use the Trocar System with appropriate locking Drill Guide to drill the pilot holes, utilizing the proper technique as previously described. The Drill Guide is threaded to the plate to ensure perpendicular pilot hole drilling.

(Reference General STG pg 22, Step 7)

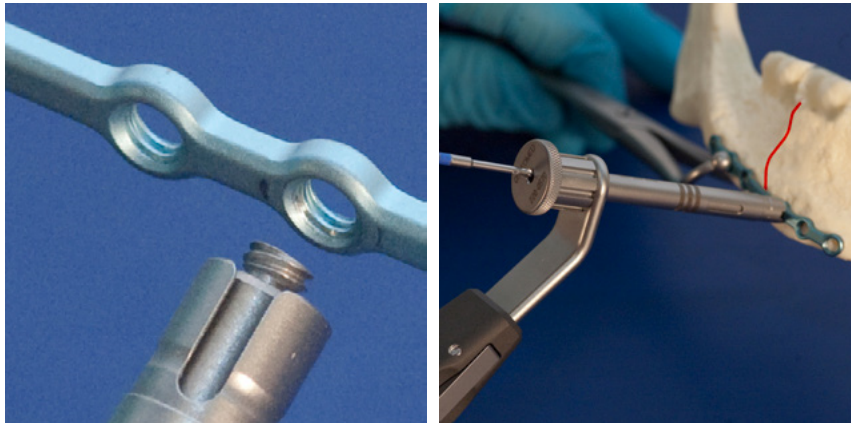




5. Use the Depth Gauge (220-0572) to determine the appropriate screw length.

(Reference transbuccal approach for Cheek Retractor pg 30, step 7)

Locking Screws

Surgical Technique Guide

6. Drive the first screw into proximal pilot hole until the screw is seated into the plate. Drive the second screw into the pilot hole on the opposing side of the fracture line/osteotomy until the screw is locked into the plate. Drive the remaining screws in a distal direction from fracture line/osteotomy.

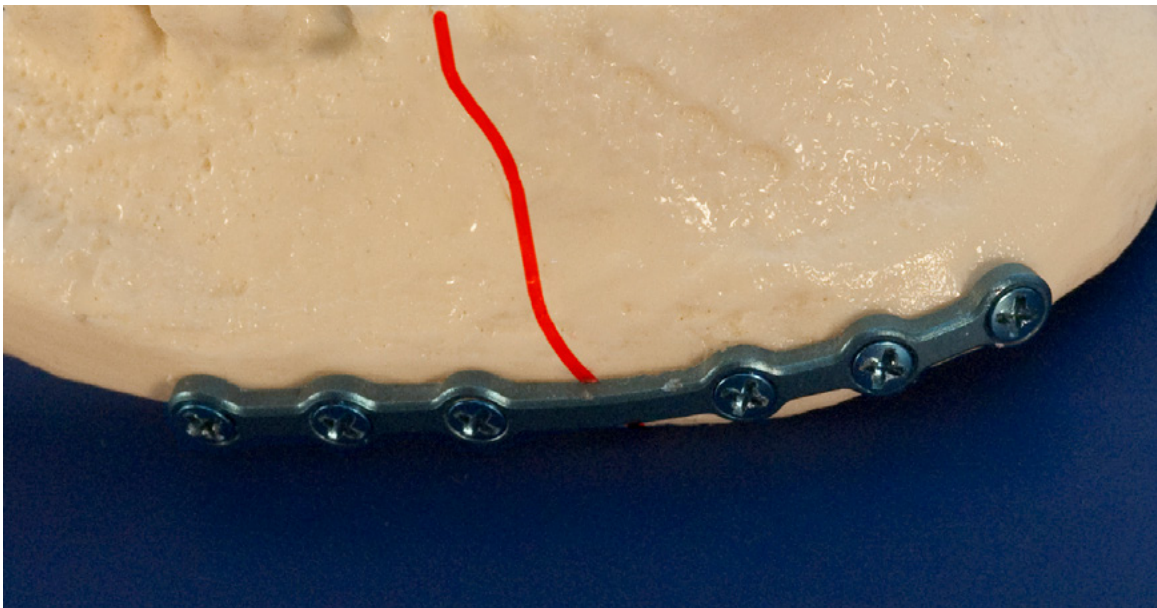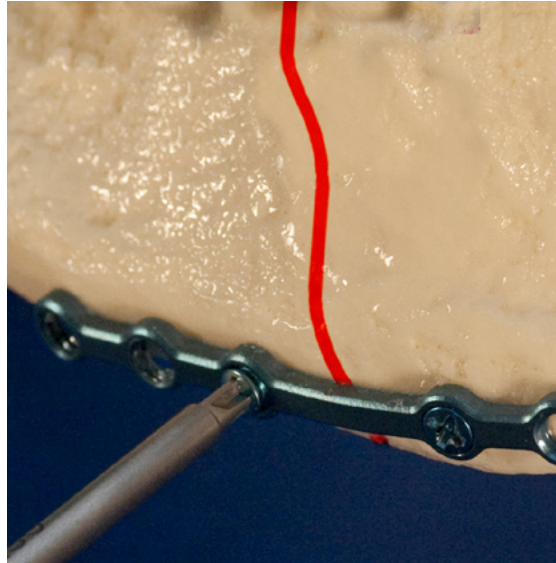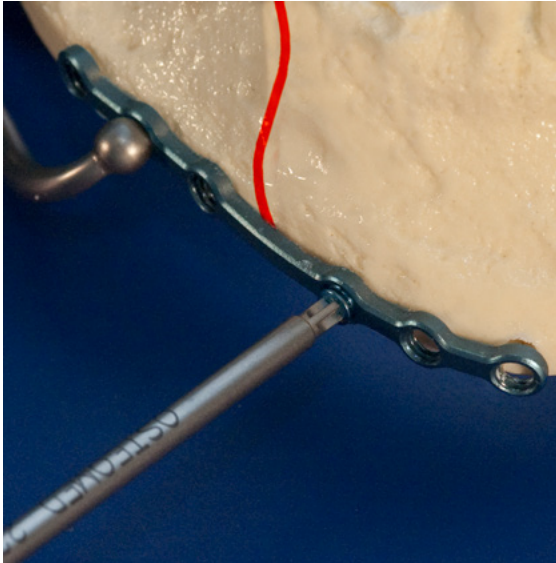

Angulated Locking Screws

Surgical Technique Guide



1. Expose the defect
2. Select and contour the appropriate template
3. Select and contour appropriate angulated locking plate, utilizing the MFX™ contouring instrumentation (**Reference General STG pg 17-20, Step 5**)

Technique Tip: locking inserts are not necessary


Note: 2.0/2.4 mm Angulated Locking Plates are rose red

Note: Angulated Locking Auto-Drive® Screws are pink

Note: Angulated Locking Standard Screws are bronze

Note: Angulated Locking Safety Screws are gray

4. Cut the Angulated Locking Plate if necessary.
Use the Round Reconstruction Plate Cutter Assembly (220-0584, 220-0585).
(**Reference General STG pg 15-16, Step 4**)
5. Select first hole to be drilled and place Angulated Locking Drill Guide flush and perpendicular to the plate

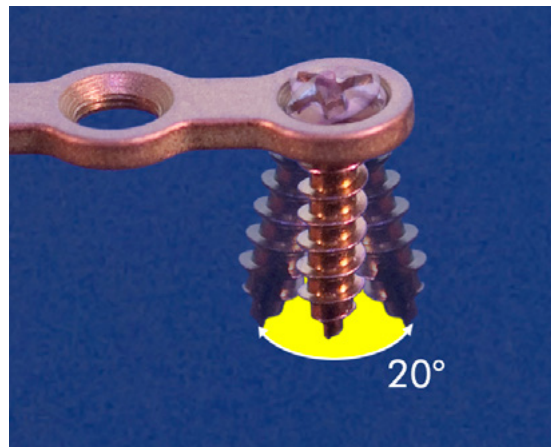

Technique Tip: Select 2.0 mm Angulated Locking Drill Guide (220-0452) or 2.4 mm Angulated Locking Drill Guide (220-0453) depending on screw/plate selection

Technique Tip: Ensure Drill Guide cone is placed perpendicular and flush to the plate. Drill Guide will provide up to 10° degrees of angulation in any direction.

Technique Tip: If transbuccal approach is needed ensure only Neutral Drill Guide is used

Angulated Locking Screws Surgical Technique Guide

6. Drill the first hole proximal to the fracture line




7. Determine screw length

Use the Depth Gauge (220-0572) to determine the appropriate screw length.

(Reference General STG pg 23, Step 8)

8. Select the first screw

Technique Tip: 2.0 mm Angulated Locking, Angulated Locking Auto-Drive® and Safety Angulated Locking screws require 1.6 mm driver. The 1.6 mm driver has a white band.

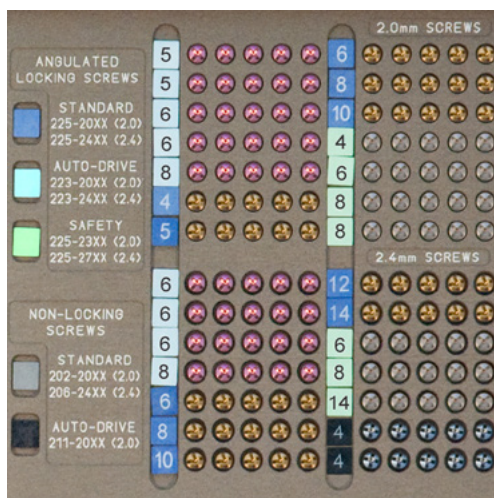


Angulated Locking Screws


Surgical Technique Guide



- 9.** Drive the screw
Insert the first screw into proximal pilot hole until the screw is seated into the plate.

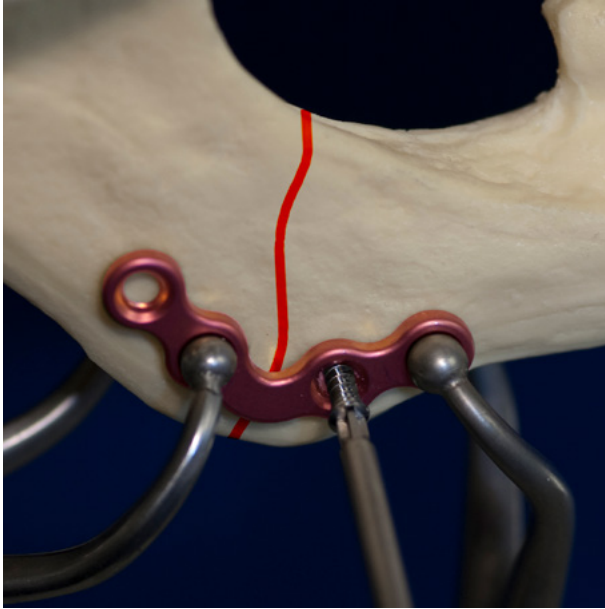

- 10.** Insert remaining screws
Drive the second screw into the pilot hole on the opposing side of the fracture line/osteotomy until the screw is locked into the plate. Repeat screw insertion steps until all screws are placed in final fixation is achieved

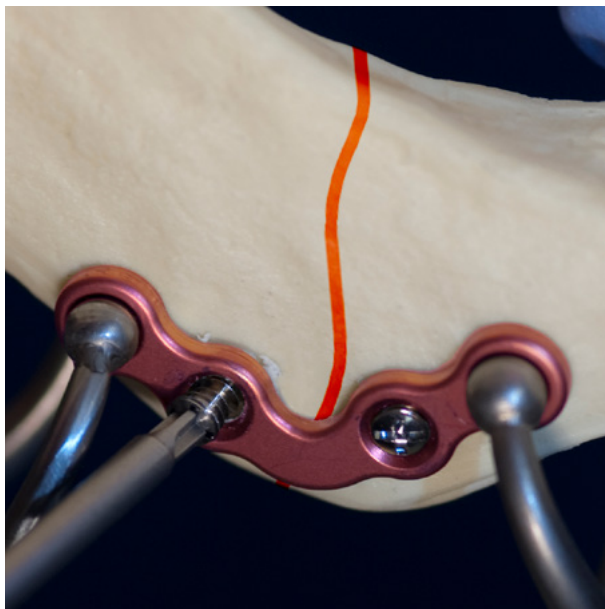

Angulated Locking Screws

Surgical Technique Guide


Note: *When placing additional screws, ensure that adjacent screw placements do not interfere with one another.*

11. Close per standard practice


Reduction Plates and Forceps Surgical Technique Guide

For mandibular symphysis, parasymphysis and body fractures only

1. Expose the fracture
2. Select the appropriate plate

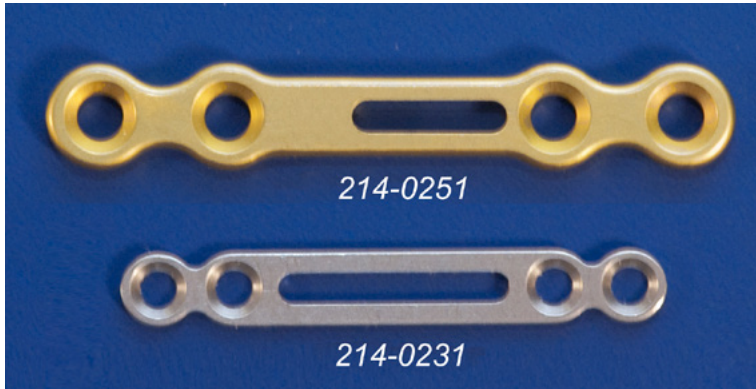

3. Select the corresponding template

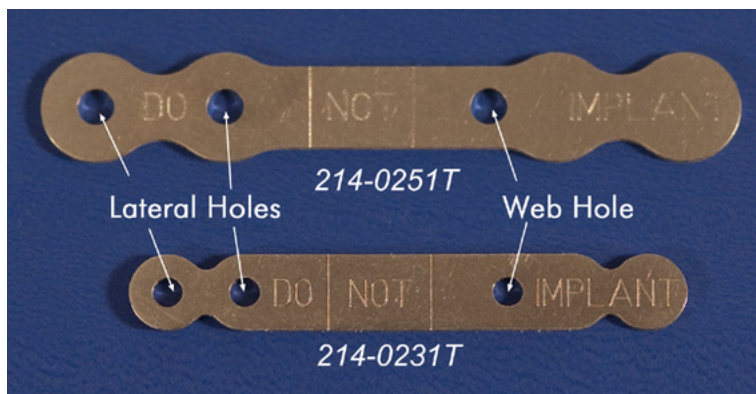


4. Drill forceps insertion holes
 - Place the template across the fracture; ensure that the fracture lines fall entirely within the limit lines on the template.
 - Drill the first insertion hole through the template's web hole and the second hole through one of the lateral holes.


Reduction Plates and Forceps

Surgical Technique Guide

Note: The etched lines present on the template indicate the maximum fracture displacement that can be reduced using this technique. If the fracture spans past the limit lines, conventional reduction methods must be used.

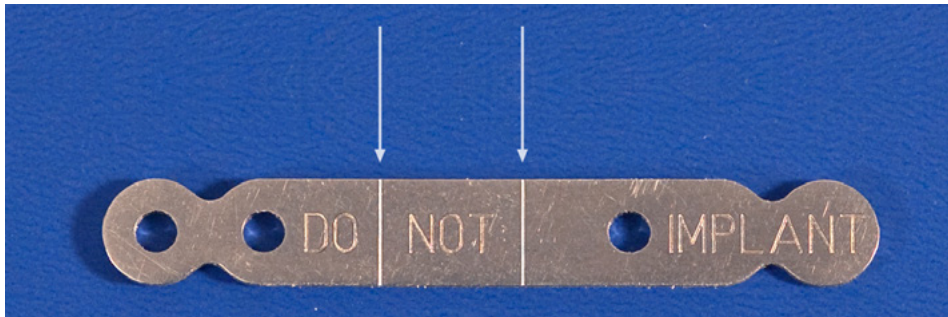

Note: The lateral hole will align with the screw hole on the plate and serve as a pilot hole during screw placement. The other hole will fall within the slot on the plate.

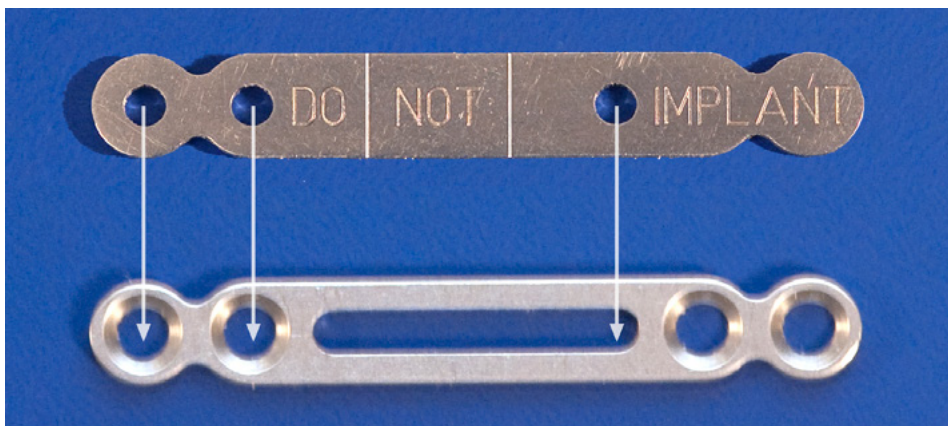

Note: The holes will be drilled with the pilot drill corresponding to the screws used with the plate.

Note: Ensure that the template does not move location during the drilling of the holes. Ensure holes are drilled as close to parallel as possible.



Reduction Plates and Forceps Surgical Technique Guide

5. Remove the template

Technique Tip : If little or no contouring of the plate is required go to Step 6, if significant contouring is required continue to Step 5.1

5.1. Reduce the fracture

Insert the forceps into the pre-drilled holes on the fragments and reduce the fracture.

5.2. Contour the template by placing above or below the forceps

Note: *Plate must be flush with the bone to provide adequate contouring*

5.3. Contour the plate to match the template.

Technique Tip: Utilize one of the OsteoMed plate contouring options.
(Reference General STG pg 17-20, Step 5)

Note: *Multiple bending of the plate may weaken the plate and could result in implant fracture and failure.*

5.4. Remove the forceps from the bone

6. Insert the forceps' tips into the plate until there is enough engagement between the tip of the forceps and the slot of the plate to adequately retain the plate.




6.1 One tip of the forceps must be placed in the slot and the second tip must be placed in one of the screw holes on the opposing side of the fracture.

Note: *Do not use excessive force when lodging the tip into the slot*


Reduction Plates and Forceps

Surgical Technique Guide

7. Engage the reduction forceps into the mandible fracture fragments, reduce the fracture.
Note: *Forceps will continue to provide reduction/compression during screw placement.*
8. Push the plate flush with the bone and drill the first screw hole using the appropriate drill.

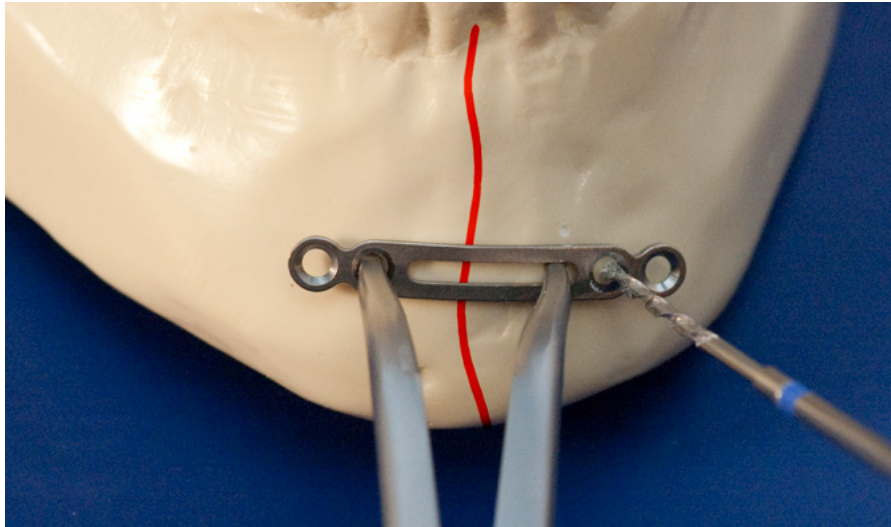

9. Use Depth Gauge to measure the screw length required.
(Reference General Section pg 23, Step 8)
10. Insert the required screw and drive until fully seated.

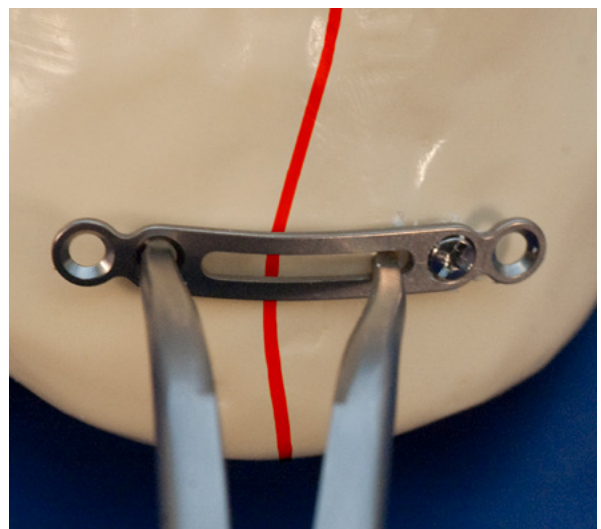


Reduction Plates and Forceps Surgical Technique Guide



11. Place opposite side lateral screw to maintain reduction.

12. Remove forceps and place remaining screws.


13. Close per standard practice.

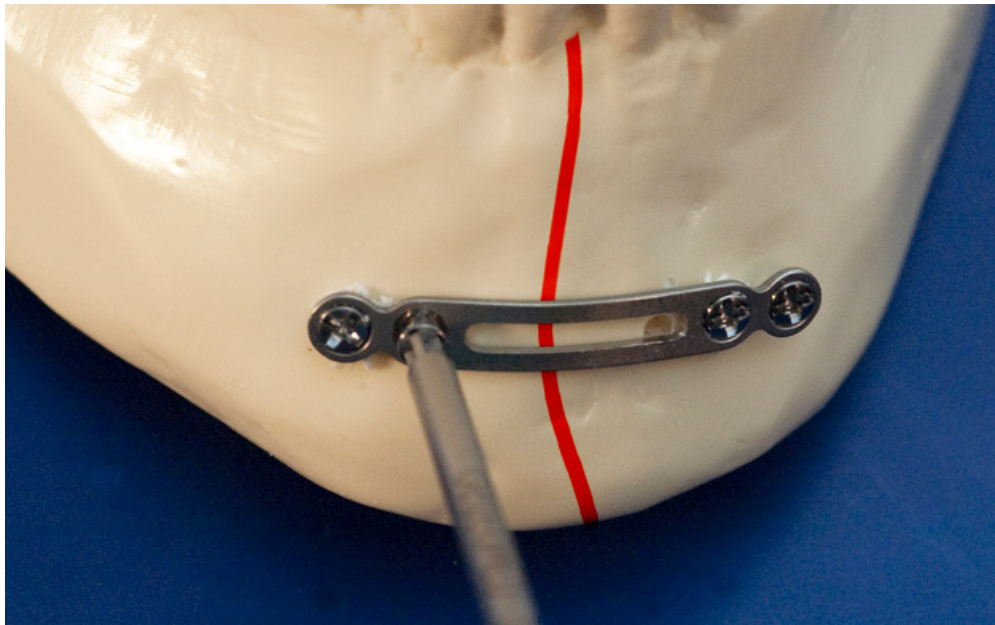


Bone Resection (Ablative) Technique

Surgical Technique Guide

1. Expose the defect and select and contour appropriate template.
2. Cut the reconstruction plate using the Round Reconstruction Plate Cutter Assembly. (220-0584, 220-0585)

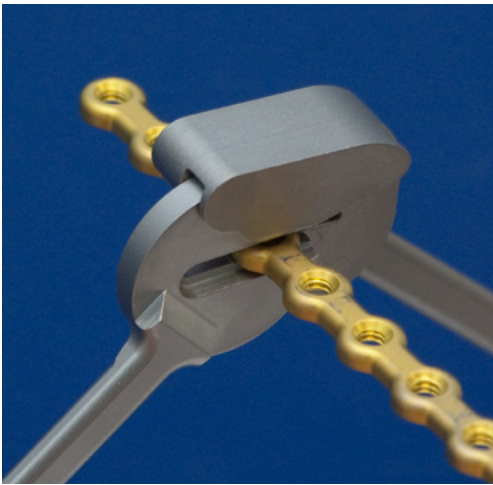

3. Place bending inserts into the reconstruction plate and contour using reconstruction plate contouring options to match the template. After contouring is complete remove Bending Inserts.

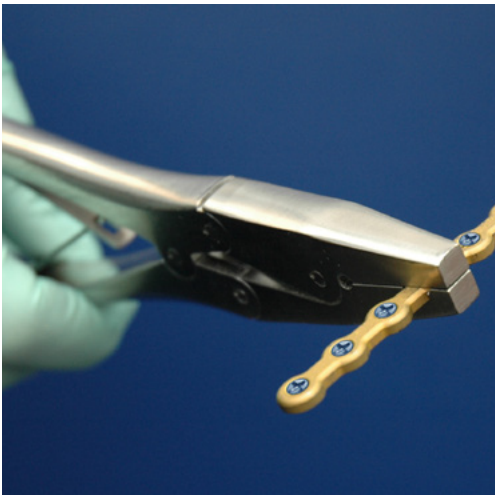


(Reference General STG pg 17-20, Step 5)

4. Drill pilot hole

Bone Resection (Ablative) Technique

Surgical Technique Guide

5. Use the Depth Gauge (220-0572) to determine the appropriate screw length and drive screws

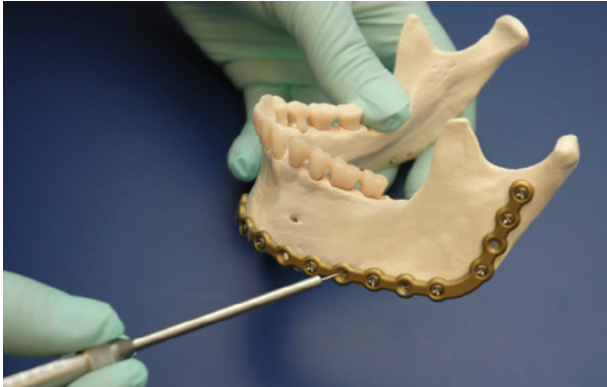

6. Once plate is in place, remove screws and plate. After mandibular resection, replace plate and screws.

Note: Note screw location upon removal. Screws need to be returned to previous locations.

7. Resect the mandible and replace implant



8. Resection is complete (If necessary apply a bone graft)



Lag Screw Device

Surgical Technique Guide

Indications

Indicated for symphyseal fracture fixation in the mandible

Lag Screw Device Instrumentation

- Modular Handle, Cannula: 220-0578
- MFX™ Lag Screw Cannula: 220-0724
- Sliding Depth Gauge, Lag Screw Depth Gauge: 220-0538
 - Length indicator slides along Lag Screw Cannula to determine screw length.
 - Utilized to predict exit point of drill on far cortex
- 2.4 mm Lag Screws 20–38 mm: 306-24XX
- Housed in the MFX™ Sterilization Tray: 220-0720


Lag Screw Device Surgical Technique Guide

1. Reduce the fracture


Instrument Tip: The Bone Reduction Forceps (220-0231) may be necessary to reduce the fracture fragments. Pre-drill insertion holes for the forceps on either side of the fracture.

2. Assemble the Lag Screw Device

- Insert Lag Screw Cannula(220-0724) into Modular Cannula Handle (220-0578)
- Glide the Lag Screw Depth Gauge (220-0538) over the cannula

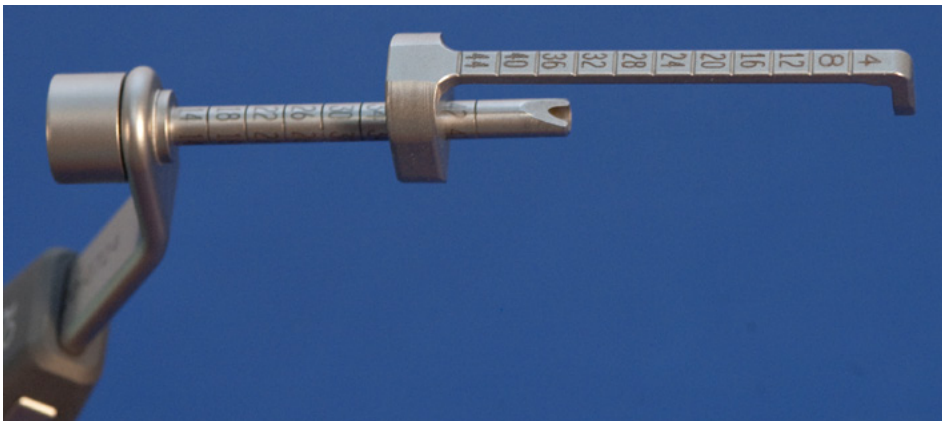

3. Drill the pilot hole

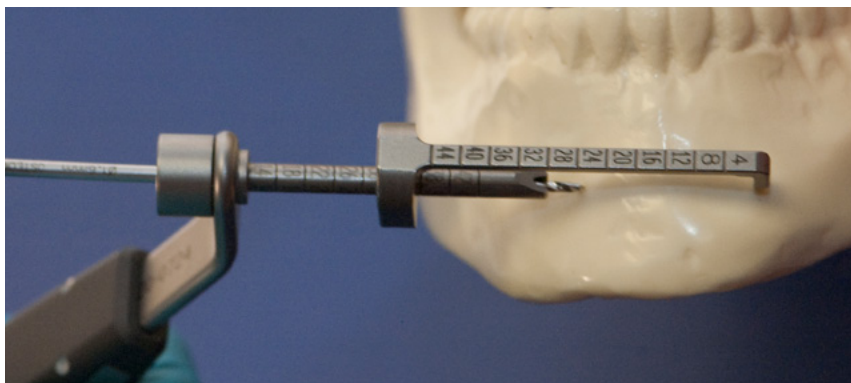

Instrumentation Tip: Once device has been completely assembled, place lag screw device onto the mandible. Ensure device is placed perpendicular to the fracture line. If desired length is not achieved, adjustment of the device can be achieved by gliding the depth gauge along the cannula. This will provide the desired exit point of the screw.

Caution: When utilizing this device make certain the device is placed in a fashion that avoids the path of the tooth roots.

Lag Screw Device

Surgical Technique Guide

Instrumentation Tip: Utilize either 220-0560 or 220-0563


- 220-0560, 2.0 mm Pilot Drill, Long, Mandible System, J-latch
- 220-0563, 2.0 mm Pilot Drill, Long, Mandible System, Drill bit, OsteoPower

Instrumentation Tip: Lag Screw Device must contact the bone at both the cannula and lag screw pointer tips. Drill entry point should be located approximately 12–15 mm from the fracture line. Drill through the first cortex past the fracture line and through the distal cortex, utilizing the long mandibular drill bit and the MFX™ Lag Screw Cannula.

4. Determine length of lag screw

Option 1: Sliding Depth Gauge (220-0538) on Lag Screw Cannula

Ensure that the tip of sliding length gauge is fully contacting symphysis before reading

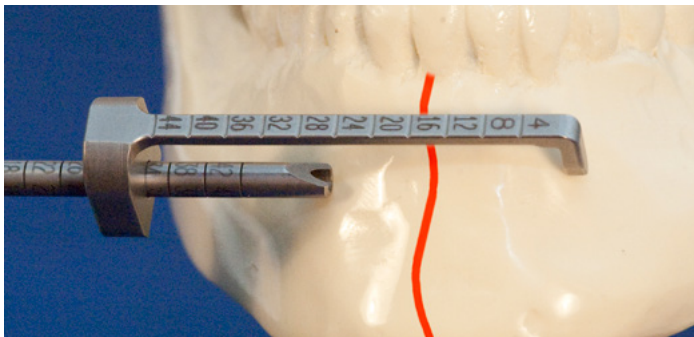

Option 2: Depth Gauge (220-0712)

5. Insert screw

Note: To prevent fragment rotation a secondary screw is recommended.

To place additional screw repeat steps 3 – 5, placing additional screw parallel to the initial screw.

MMF Screw System



System Components

- 2.0 mm MMF Screw (209-20XX) AutoDrive®
Lengths: 8 mm, 11 mm, 14 mm
- 2.4 mm MMF Screw (207-24XX) Quick-Fix™
Lengths: 12 mm, 14 mm, 16 mm, 18 mm
- Check on driver stem
- 2.0 mm/2.4 mm Screwdriver Shaft, Manual (220-0021)
- TaperLock™ Screwdriver Body (220-0019)
- 2.0 mm Pilot Drill, J-latch (220-0005)
- 2.0 mm Pilot Drill, Manual (220-0008)
- 24 ga. Stainless Steel Wire (207-0120)

Features

- Streamlined system eliminating the need for arch bars, dramatically reducing application time for MMF
- Smooth screw posts minimizing soft tissue irritation
- Reduced intraoral hardware increasing patient comfort
- Wire passing holes providing secure ligature wire engagement
- TaperLock™ Cruciform Drive carries screws to the site with positive engagement

Indications

- Orthognathic procedures
- Endentulous or partially endentulous
- Simple fracture patterns

Contraindications

- Comminuted fractures
- Unstable fracture segments
- Pediatric injuries


MMF Screw System

Surgical Technique Guide

1. Plan

Select the number and desired position of screw placement. This will be based on fracture pattern and anatomic location.

2. Locate maxillary tooth buds

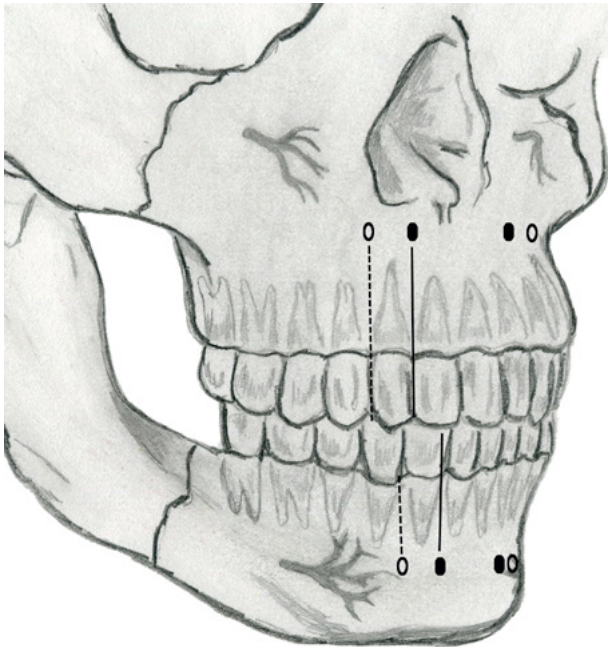

Technique Tips: Pay special attention to the canine tooth root, longest of all roots. Avoid all existing dentition while locating and avoiding the infraorbital and mental nerves. All dentition should be authenticated through the usage of the appropriate radiograph: panoramic X-ray.

Note: Maxillary MMF screws should be placed 5 mm superiorly to the tooth root.

3. Select appropriate MMF Screws

Based on fracture pattern, patient anatomy and anticipated patient compliance select either 2.0 mm (209-20XX) or 2.4 mm (207-24XX)

MMF Screw System

Surgical Technique Guide



4. Insert first MMF screw into the maxilla



Technique Tip: Screws may be placed through the mucosa without an incision. Insert either the 2.0 mm or 2.4 mm MMF screw into the maxilla, paying close attention not to apply unnecessary pressure to the mucosa. Drive the screw to the appropriate depth, leaving the wire-passing hole exposed. Do not over torque or bottom out the screw.

Note: *When utilizing 2.4 mm screws or when dense bone is encountered create a pilot hole*

Note: *2.0 mm MMF Auto-Drive® screws are self-drilling*

5. Insert MMF screw into the mandible

Start by identifying the necessary anatomic landmarks: tooth roots and mental nerve.

6. For secondary screw ensure placement in the mandible is 5 mm inferior and medial or lateral to the canine tooth roots.

7. Insert additional screws


Follow the previously outlined procedure to place the remaining screws.



Note: A minimum of 3 pairs of MMF screws are recommended to ensure adequate stability. A pair consists of one screw in the mandible and an opposing screw in the maxilla.

8. Insert Wire

Insert 24 gauge stainless steel wire (207-0120) through exposed wire passing holes into maxillary and opposing mandibular MMF screw heads in a vertical and “X” pattern. Tighten only enough to provide provisional fixation.

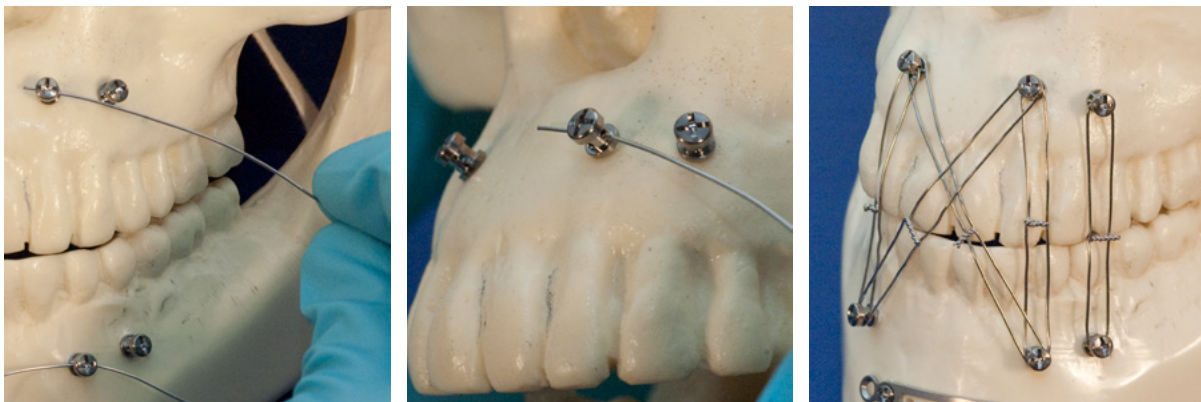

Note: Alternatively, wires may be wound around grooves in screw heads.

MMF Screw System

Surgical Technique Guide



9. Establish occlusion

Note: *Prior to securing all MMF screws, occlusion must be established*

10. Tighten wires fully



www.acumed.net

Acumed USA Campus
5885 NE Cornelius Pass Road
Hillsboro, OR 97124
+1.888.627.9957

OsteoMed USA Campus
3885 Arapaho Road
Addison, TX 75001
+1.800.456.7779

Acumed Iberica Campus
C. Proción, 1
Edificio Oficor
28023 Madrid, Spain
+34.913.51.63.57

CMF00-00-A | Effective: 2024/06 | © 2024 Acumed® LLC

These materials contain information about products that may or may not be available in any particular country or may be available under different trademarks in different countries. The products may be approved or cleared by governmental regulatory organizations for sale or use with different indications or restrictions in different countries. Products may not be approved for use in all countries. Nothing contained in these materials should be construed as a promotion or solicitation for any product or for the use of any product in a particular way that is not authorized under the laws and regulations of the country where the reader is located. Nothing in these materials should be construed as a representation or warranty as to the efficacy or quality of any product, nor the appropriateness of any product to treat any specific condition. Physicians may direct questions about the availability and use of the products described in these materials to their authorized Acumed distributor. Specific questions patients may have about the use of the products described in these materials or the appropriateness for their own conditions should be directed to their own physician.

Refer to the provided instructions for use for the complete indications, contraindications, warnings, and instructions for use.

OsteoMed® LLC is a wholly owned subsidiary of Acumed LLC.
OsteoMed and MFx™ are registered trademarks of OsteoMed LLC.
Acumed® is a registered trademark of Acumed LLC.